Microtech 648ISO World

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The new Microtech 648ISO is powered by the new Software Elfolab capable to satisfy all operator's needs, thanks to its remarkable versatility. By connecting the PC to a Host Computer, Elfolab can easily import and export databases. Elfolab is extremely user friendly. Making data storage, customized analysis set up, 4 different reports, constant instrument monitoring has never been so easy for the operator. Thanks to the combination of the new Microtech 648ISO and new Elfolab, Interlab has made a giant step ahead in the Electrophoresis world.

www.interlab-srl.com

Excellence in Automated Electrophoresis

Interlab S.r.l.
Via Rina Monti, 26 - 00155 Roma (Italy)
Tel. +39 - 06 - 227.54.350 Fax +39 - 06 - 227.54.534
E-mail: info@interlab-srl.com <http://www.interlab-srl.com>

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Let the Microtech 648ISO run and you will see that it has no limits. The new Microtech 648ISO is the only instrument in the world performing in full automation also the most complex routine electrophoresis analysis such as: Isoenzymes, Unconcentrated Urine Proteins, and CSF. Total automation is Interlab's goal for the new Microtech 648ISO. We made it and we are proud of it!

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The new Microtech 648ISO can work both in micro or semimicro technique with its two applicators, 8 (micro) or 4 (semimicro) independent blades.

Up to 48 (micro) or 24 (semimicro) samples can be loaded at once and kept cold due to the Peltier effect.

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
The monobloc migration chamber with the two electrodes (carbonium) allows a perfect current distribution over the entire strip surface.

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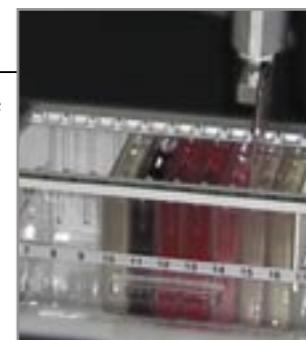The new Multi Reagent Independent Tanks is a reagent loading system that allows the operator a quick replacement of the reagents. In addition, the low reagent consumption and the absence of water source (pumps, tubes etc.) make the system highly reliable.

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
The heart of the instrument is the intelligent electromagnetic head on which a built-in optical sensor ensures always high precision and reliability.

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
The revolutionary Microtech 648ISO way of reading, with its integrated 8 channels densitomer, guarantees speed, precision, and reliability. The new incubator, integrated in the system, is the key point for the automated execution of the Isoenzymes analysis.

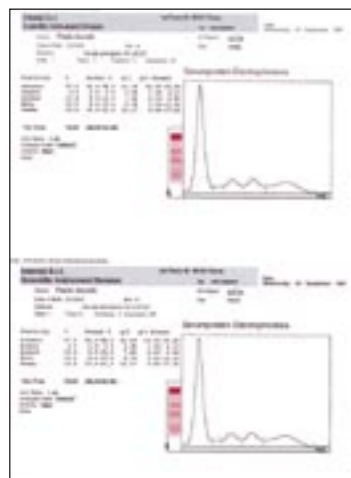Sieroproteine
Serum proteins

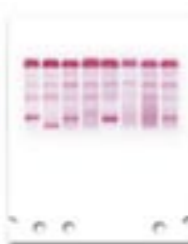
Emoglobine
Hemoglobins

Lipoproteine
Lipoproteins

Proteine Urinarie
non concentrate
Unconcentrated
Urine Proteins

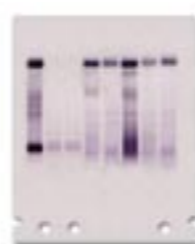LDH

ALP

CPK